

Introduction to NGS

Hubert Rehrauer




University of
Zurich UZH

ETH

Eidgenössische Technische Hochschule Zürich
Swiss Federal Institute of Technology Zurich

NGS Data Increase

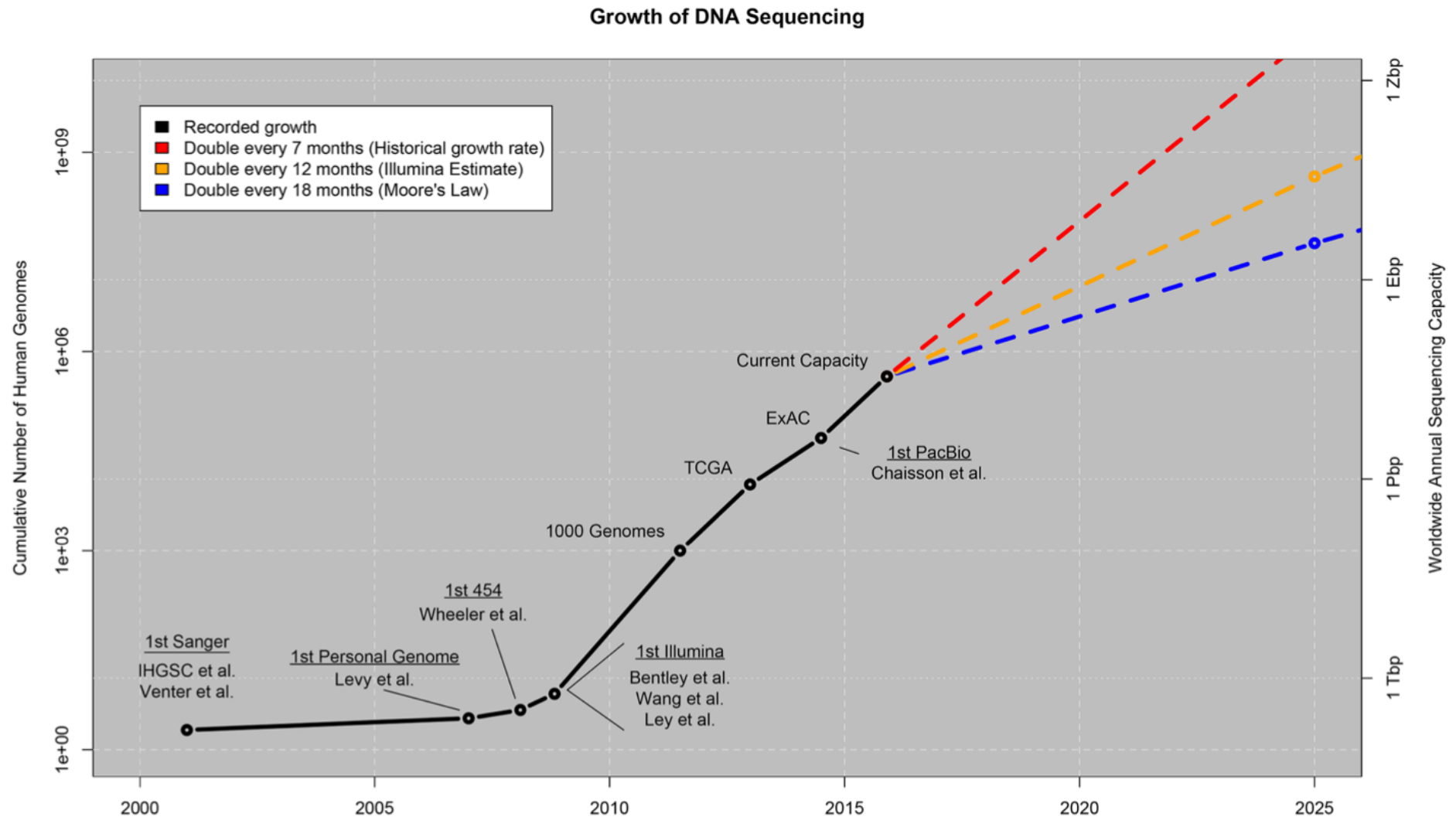

- NGS data increases faster than computer speed

Stephens ZD et al. (2015) Big Data: Astronomical or Genomical? PLoS Biol 13(7): e1002195.

Ingredients for the success

- Evolution has yielded DNA and RNA molecules for information storage and transfer. They have good properties to be read (**measured**)
- NGS technologies rely on
 - **massive parallelization**
 - measurement process is done by individual molecules (**cheap and fast**)

2nd generation sequencers

3rd generation sequencers

Illumina

Pac Bio



optical

Ion Torrent


Oxford Nanopores


non-optical

Illumina Flow cell

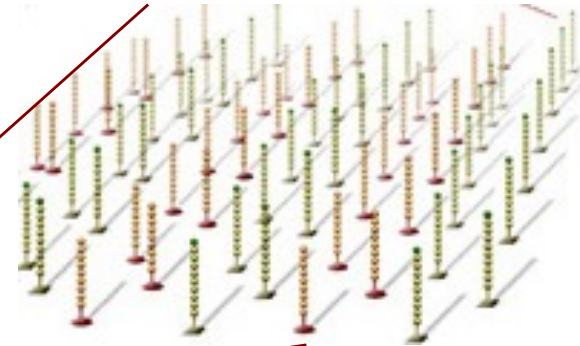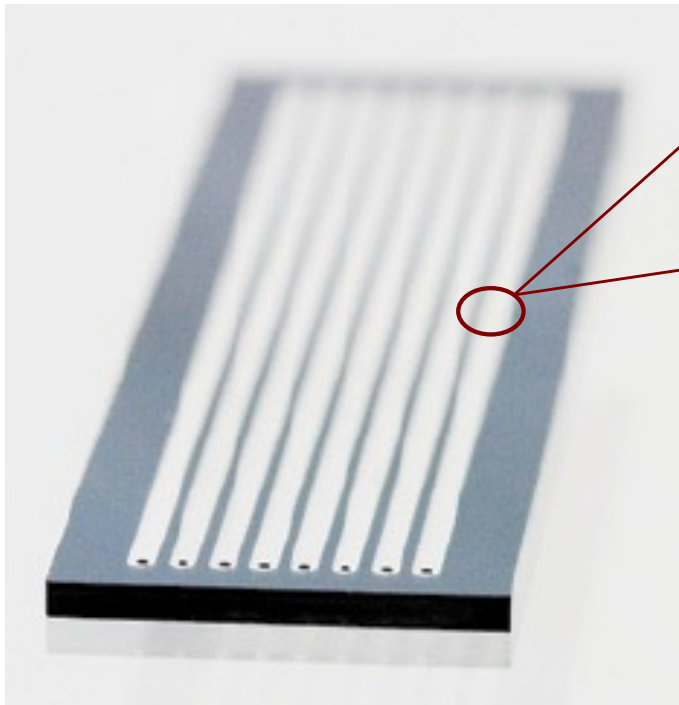



Cluster generation overview

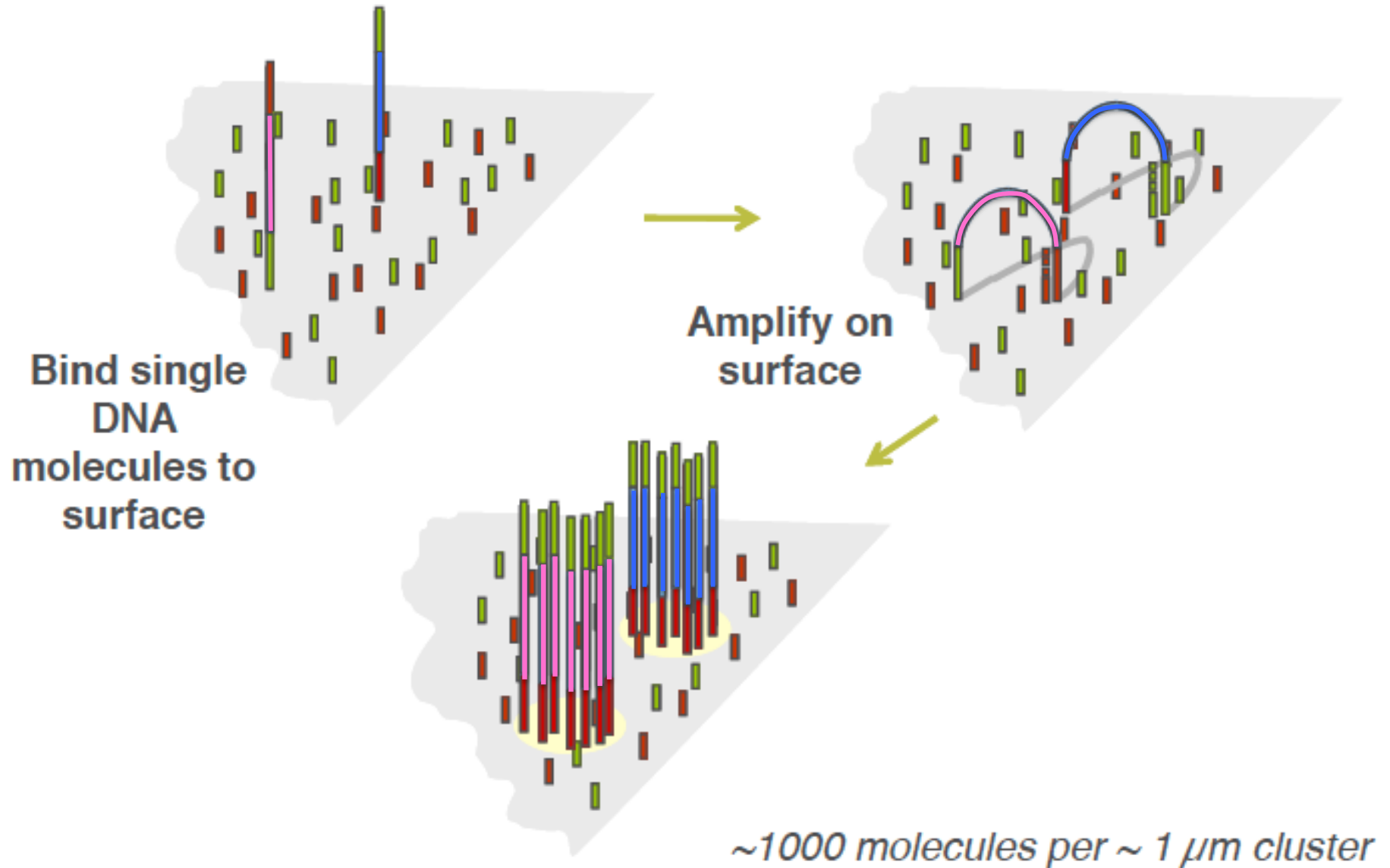


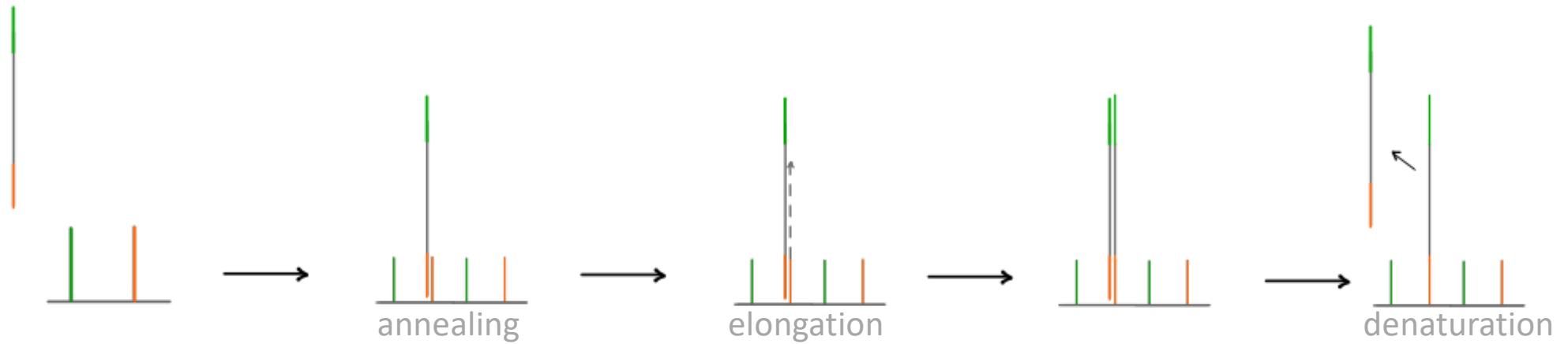


bridge amplification

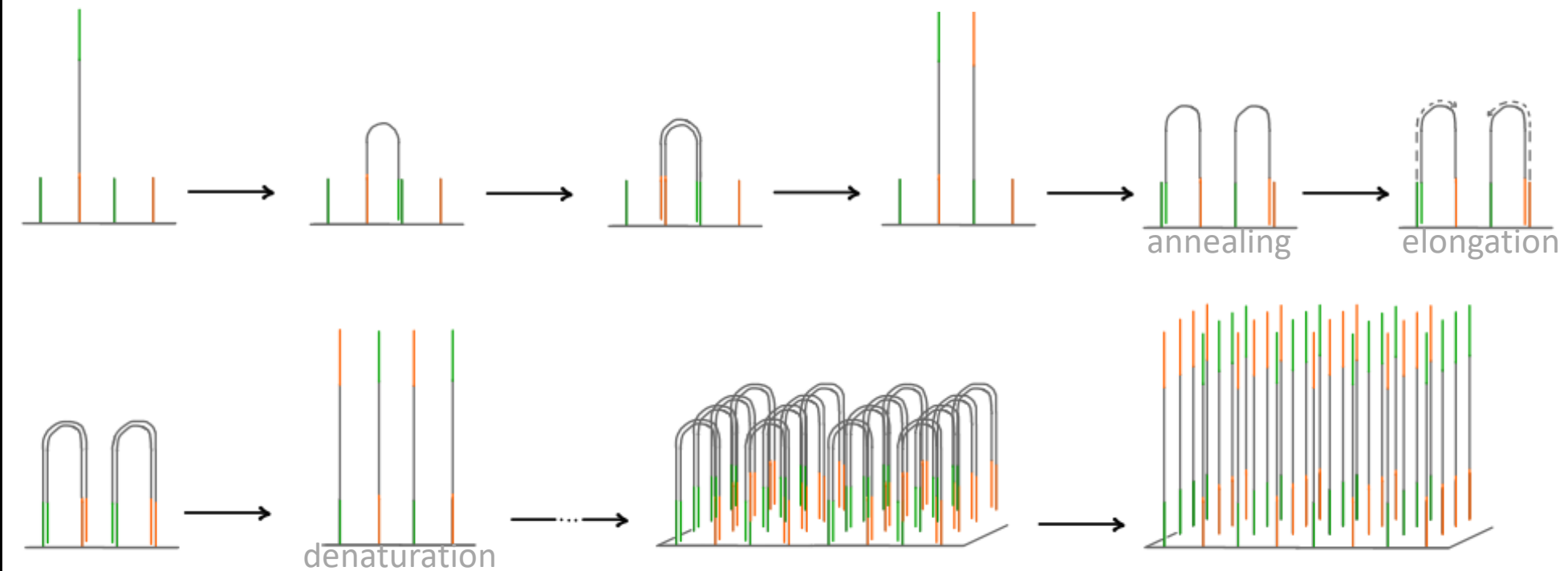


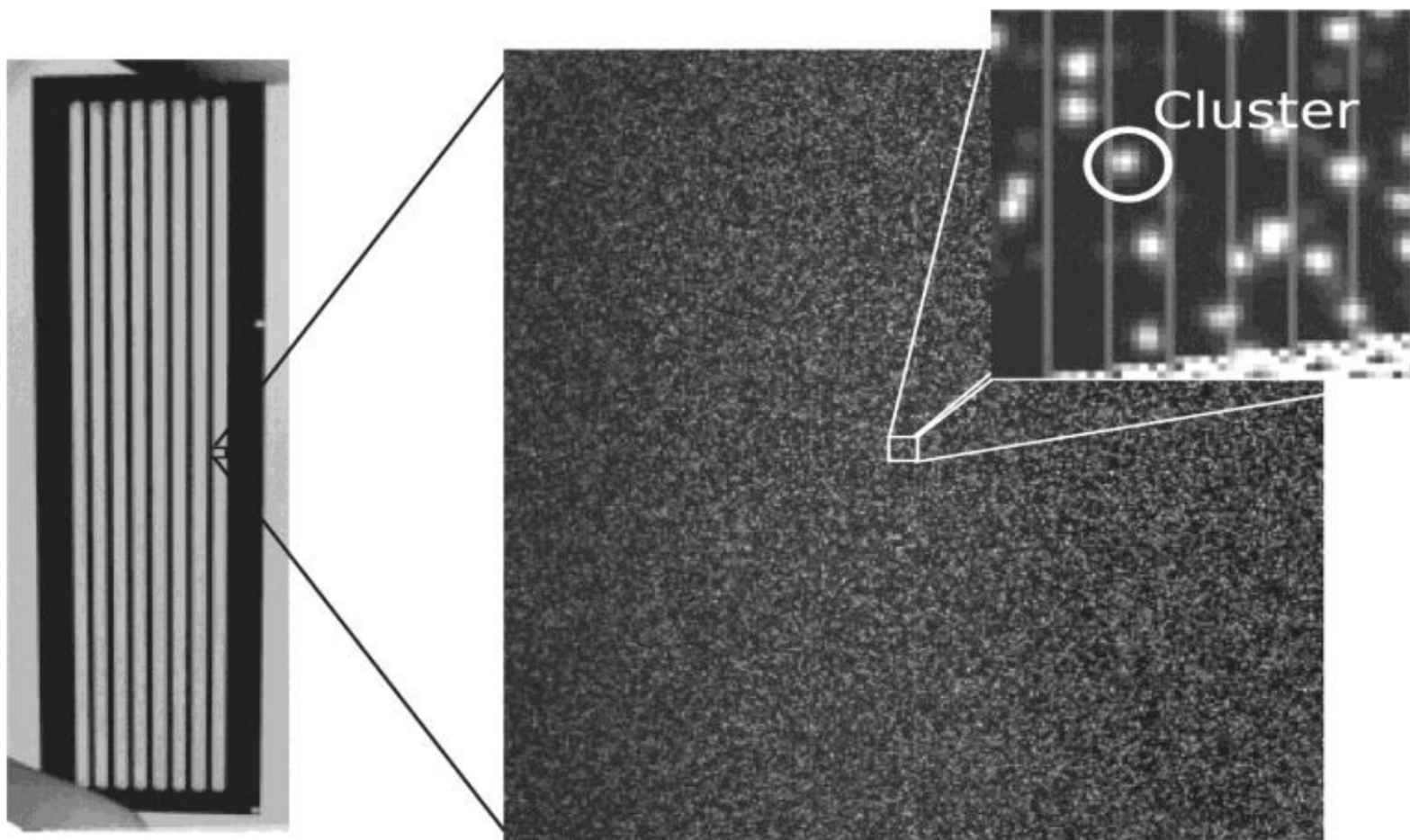


Illumina Sequencing

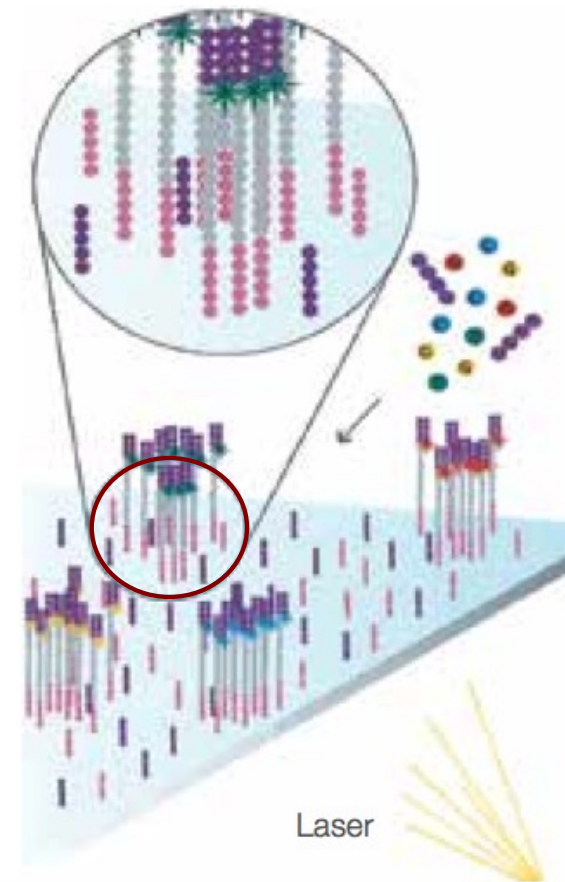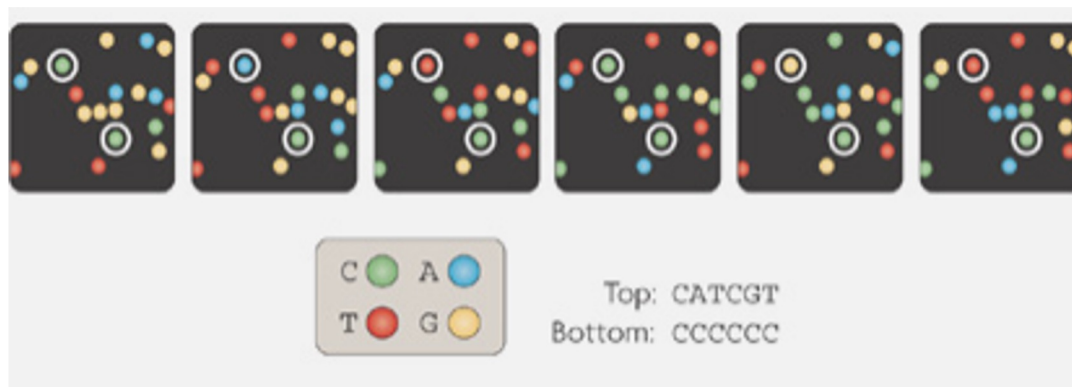


The first sequencing cycle begins by adding four labeled reversible terminators, primers, and DNA polymerase.

Color coding of bases

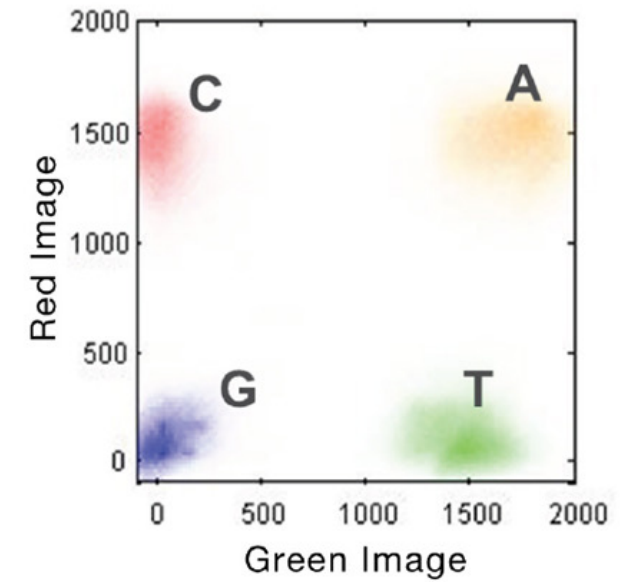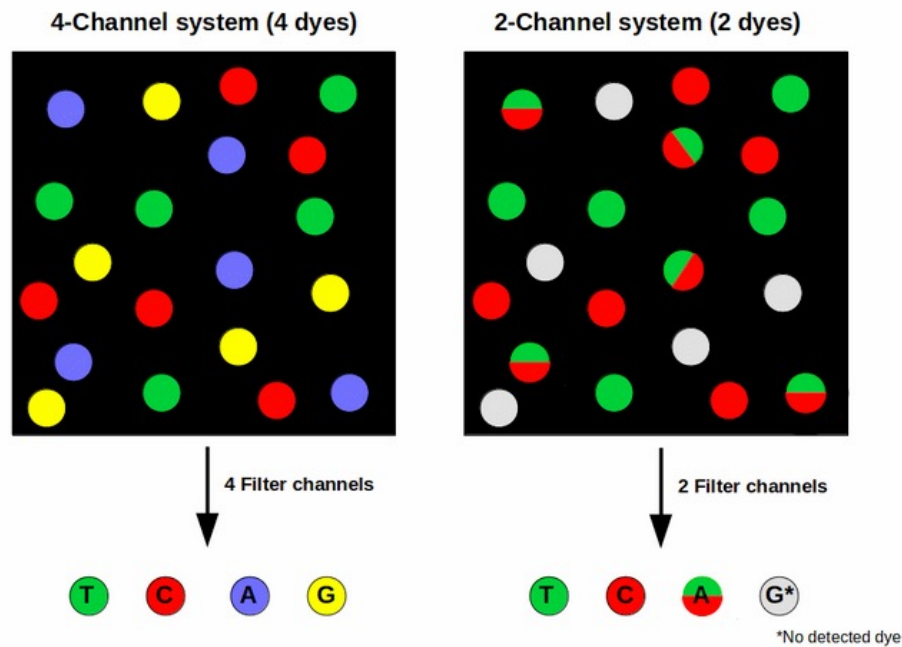

Color coding of illumina iSeq

A.

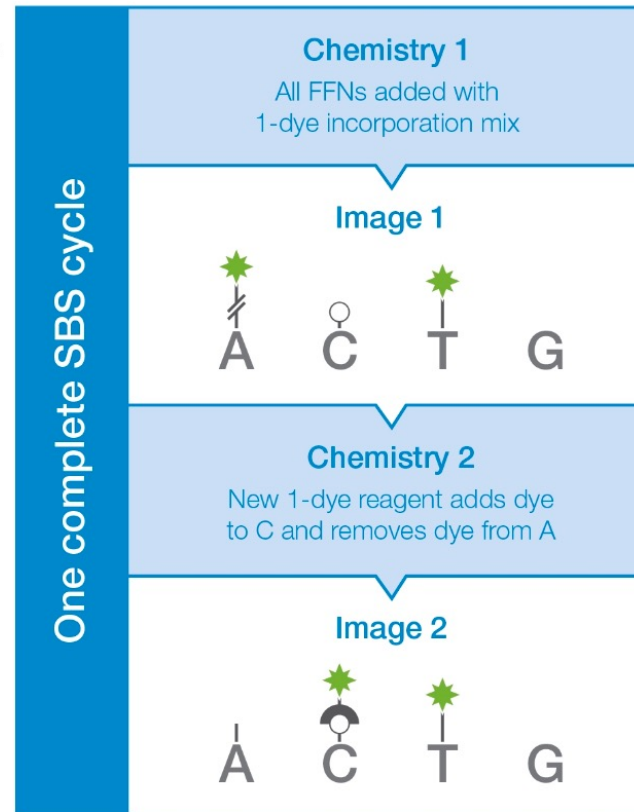

B.

Image 1	Image 2	Result
ON	OFF	A
OFF	ON	C
ON	ON	T
OFF	OFF	G

Phred scores measure base call accuracy

- P
 - error probability of a given base call
- Q
 - $-10\log_{10}P$
- Assign to each base
- Range from 0-41



Phred Quality Score	Probability of Incorrect Base Call	Base Call Accuracy
10	1 in 10	90%
20	1 in 100	99%
30	1 in 1,000	99.9%
40	1 in 10,000	99.99%
50	1 in 100,000	99.999%

Ewing B, Green P. 1998. Genome Res. 8(3):186-194.

http://en.wikipedia.org/wiki/Phred_quality_score

```
@HWI-ST1034:40:C08PJACXX:2:1101:20681:1994 1:N:0:ATCACG
CTCGNAGACTGGCAACTTGTTCTGGTTTACTGCACCTTCTTTTAAAGGCAGAAAGGC
+
CCCF#2ADHHHGHJJJIJJHIIIIJJHIJJJJJJJJJBGIIJJJJJJJJJJJJJJJJJJ
```

Phred scores can be ASCII encoded

- Add an offset and convert the sum to ASCII
- Current format
 - **Illumina 1.9 (i.e. Sanger format)**
 - Phred scoring: 0-41;
 - Offset: 33
 - $41+33=74$ (J)
 - All current sequencers

Dec	Hx	Oct	Char	Dec	Hx	Oct	Html	Chr	Dec	Hx	Oct	Html	Chr	Dec	Hx	Oct	Html	Chr
0	0	000	NUL (null)	32	20	040	 	Space	64	40	100	@	@	96	60	140	`	`
1	1	001	SOH (start of heading)	33	21	041	!	!	65	41	101	A	A	97	61	141	a	a
2	2	002	STX (start of text)	34	22	042	"	"	66	42	102	B	B	98	62	142	b	b
3	3	003	ETX (end of text)	35	23	043	#	#	67	43	103	C	C	99	63	143	c	c
4	4	004	EOT (end of transmission)	36	24	044	$	\$	68	44	104	D	D	100	64	144	d	d
5	5	005	ENQ (enquiry)	37	25	045	%	%	69	45	105	E	E	101	65	145	e	e
6	6	006	ACK (acknowledge)	38	26	046	&	&	70	46	106	F	F	102	66	146	f	f
7	7	007	BEL (bell)	39	27	047	'	'	71	47	107	G	G	103	67	147	g	g
8	8	010	BS (backspace)	40	28	050	((72	48	110	H	H	104	68	150	h	h
9	9	011	TAB (horizontal tab)	41	29	051))	73	49	111	I	I	105	69	151	i	i
10	A	012	LF (NL line feed, new line)	42	2A	052	*	*	74	4A	112	J	J	106	6A	152	j	j
11	B	013	VT (vertical tab)	43	2B	053	+	+	75	4B	113	K	K	107	6B	153	k	k
12	C	014	FF (NP form feed, new page)	44	2C	054	,	,	76	4C	114	L	L	108	6C	154	l	l
13	D	015	CR (carriage return)	45	2D	055	-	-	77	4D	115	M	M	109	6D	155	m	m
14	E	016	SO (shift out)	46	2E	056	.	.	78	4E	116	N	N	110	6E	156	n	n
15	F	017	SI (shift in)	47	2F	057	/	/	79	4F	117	O	O	111	6F	157	o	o
16	10	020	DLE (data link escape)	48	30	060	0	0	80	50	120	P	P	112	70	160	p	p
17	11	021	DC1 (device control 1)	49	31	061	1	1	81	51	121	Q	Q	113	71	161	q	q
18	12	022	DC2 (device control 2)	50	32	062	2	2	82	52	122	R	R	114	72	162	r	r
19	13	023	DC3 (device control 3)	51	33	063	3	3	83	53	123	S	S	115	73	163	s	s
20	14	024	DC4 (device control 4)	52	34	064	4	4	84	54	124	T	T	116	74	164	t	t
21	15	025	NAK (negative acknowledge)	53	35	065	5	5	85	55	125	U	U	117	75	165	u	u
22	16	026	SYN (synchronous idle)	54	36	066	6	6	86	56	126	V	V	118	76	166	v	v
23	17	027	ETB (end of trans. block)	55	37	067	7	7	87	57	127	W	W	119	77	167	w	w
24	18	030	CAN (cancel)	56	38	070	8	8	88	58	130	X	X	120	78	170	x	x
25	19	031	EM (end of medium)	57	39	071	9	9	89	59	131	Y	Y	121	79	171	y	y
26	1A	032	SUB (substitute)	58	3A	072	:	:	90	5A	132	Z	Z	122	7A	172	z	z
27	1B	033	ESC (escape)	59	3B	073	;	;	91	5B	133	[[123	7B	173	{	{
28	1C	034	FS (file separator)	60	3C	074	<	<	92	5C	134	\	\	124	7C	174	|	
29	1D	035	GS (group separator)	61	3D	075	=	=	93	5D	135]]	125	7D	175	}	}
30	1E	036	RS (record separator)	62	3E	076	>	>	94	5E	136	^	^	126	7E	176	~	~
31	1F	037	US (unit separator)	63	3F	077	?	?	95	5F	137	_	_	127	7F	177		DEL


Read Quality Control

- Library construction could introduce bias
 - Fragmentation, ligation, amplification
 - GC bias
 - Over-amplification
 - Contamination
- Sequencing errors
 - Chemical, optical, computational

Platform	Primary error	Error rate (%)
Illumina	Substitution	0.1
PacBio	Indel	12 (consensus: 1)
Oxford Nanopore	Indel	3 - 20

Per base sequence quality - FastQC

- Range of quality values across all bases at each position

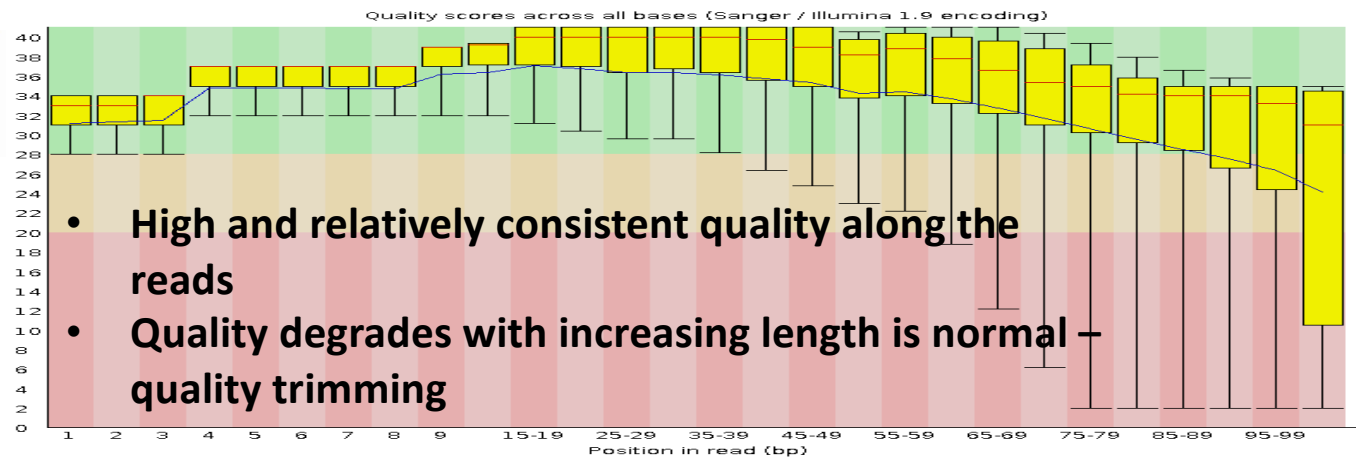


Green: >Q28, good

Orange: >Q20, reasonable

Red: <Q20, poor



Median > Q25

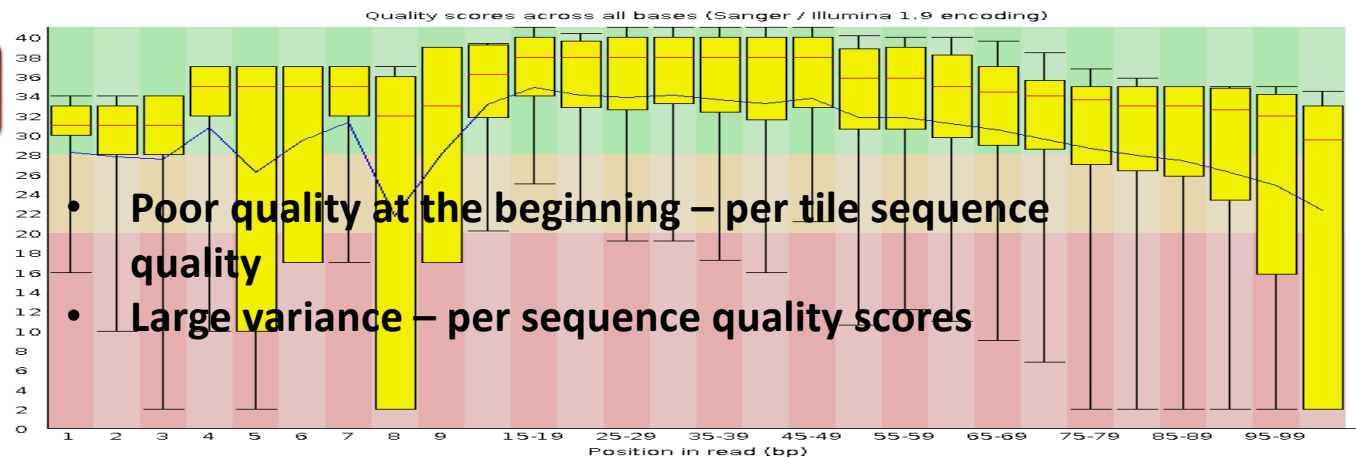

Median < Q20

Per sequence quality scores - FastQC


- Subset of sequences with universally low quality values

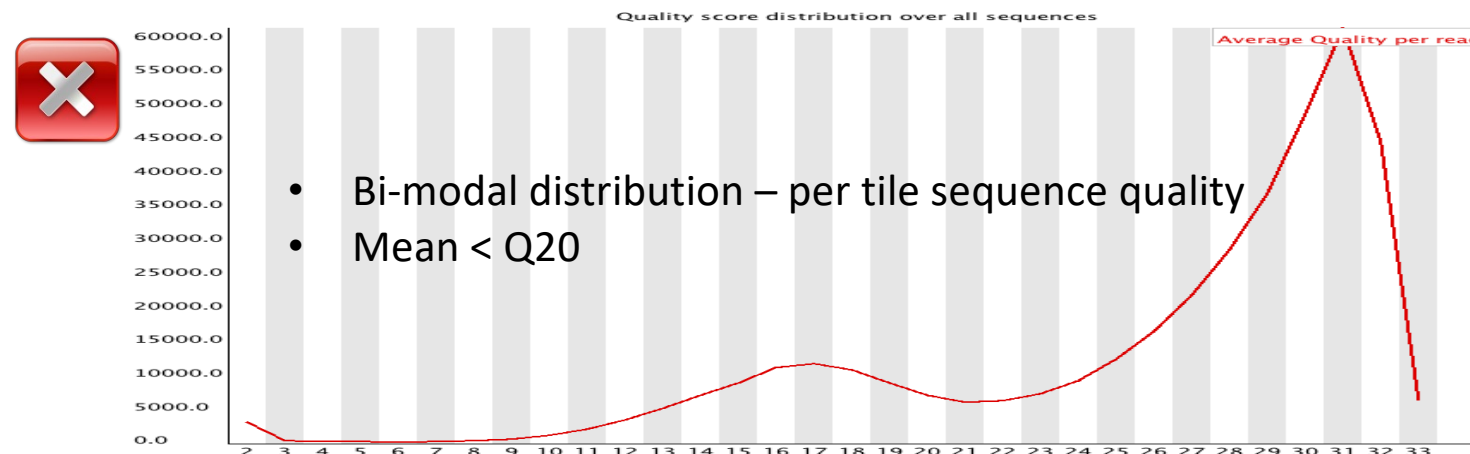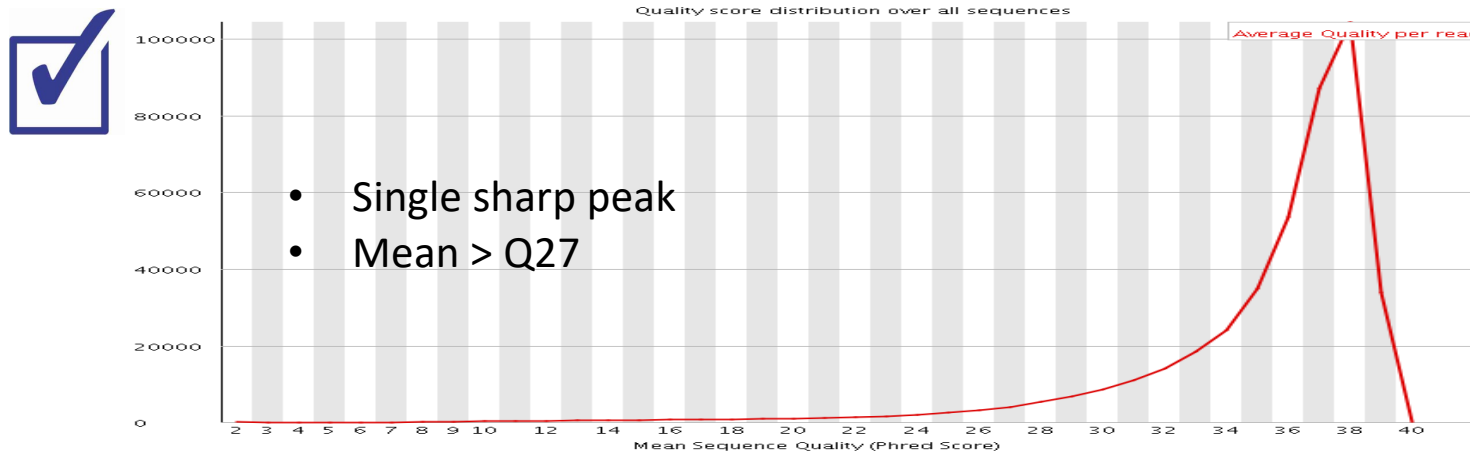

Per tile sequence quality - FastQC

- Quality scores from each tile across all bases - loss in quality associated with only one part of the flowcell

Quality per tile

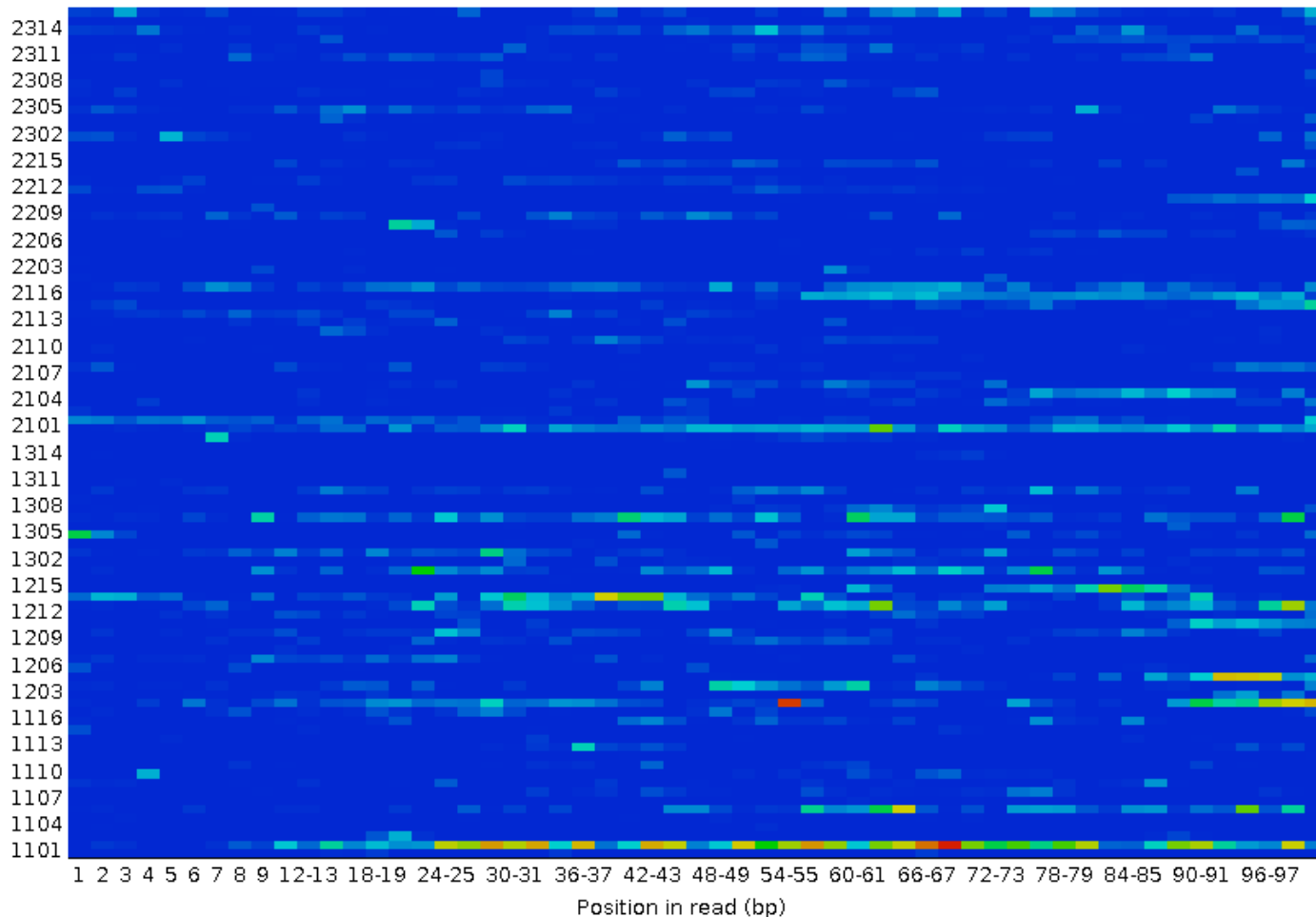

Deviation from average quality

Cold colors: \geq average




Hotter color: worse quality

Good: universal blue

Failure: $<$ average - 5

Per base sequence content - FastQC

- The portion of A, T, G, and C at each position

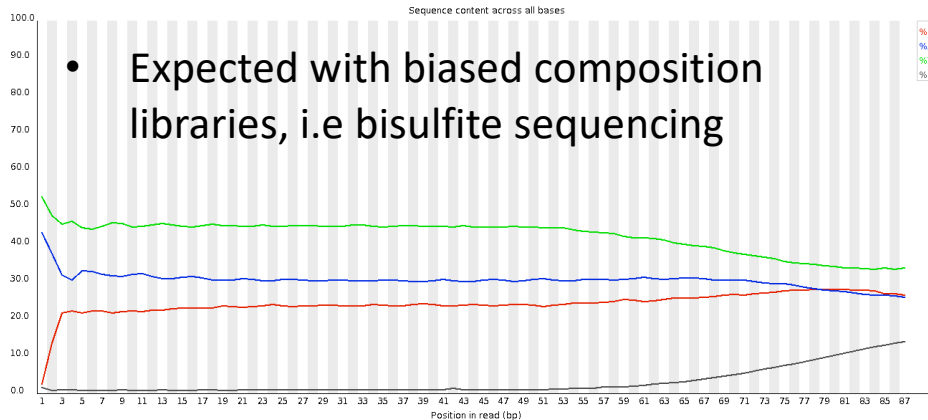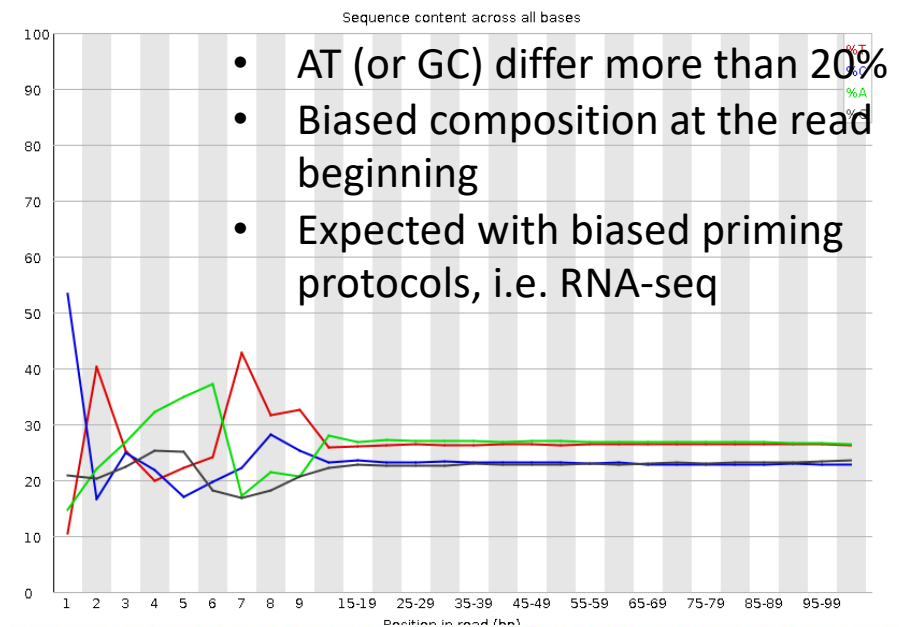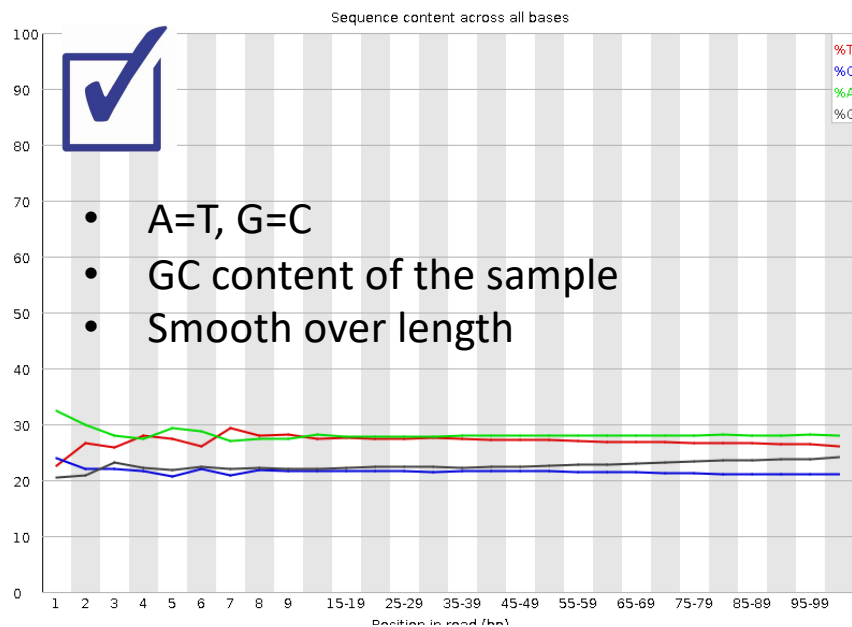


Biases in Illumina transcriptome sequencing caused by random hexamer priming

Kasper D. Hansen^{1,*}, Steven E. Brenner² and Sandrine Dudoit^{1,3}


Treatment of DNA with bisulfite converts cytosine to uracil, but leaves methylated cytosine unaffected. Therefore, DNA that has been treated with bisulfite retains only methylated cytosines.

Per sequence GC content - FastQC

- Distribution of average GC in all reads

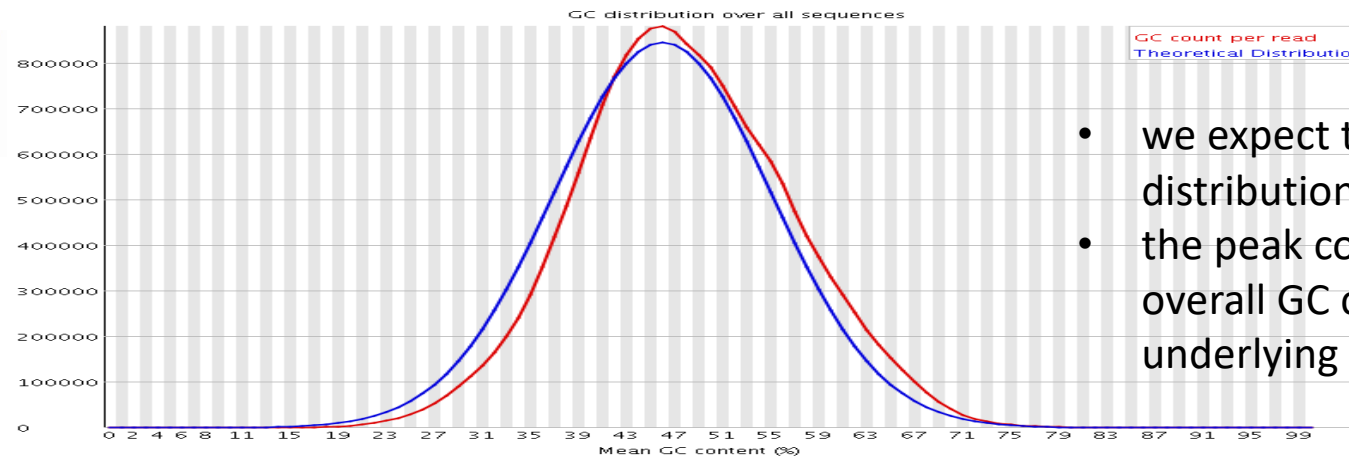


- we expect to see a roughly normal distribution of GC content
- the peak corresponds to the overall GC content of the underlying genome

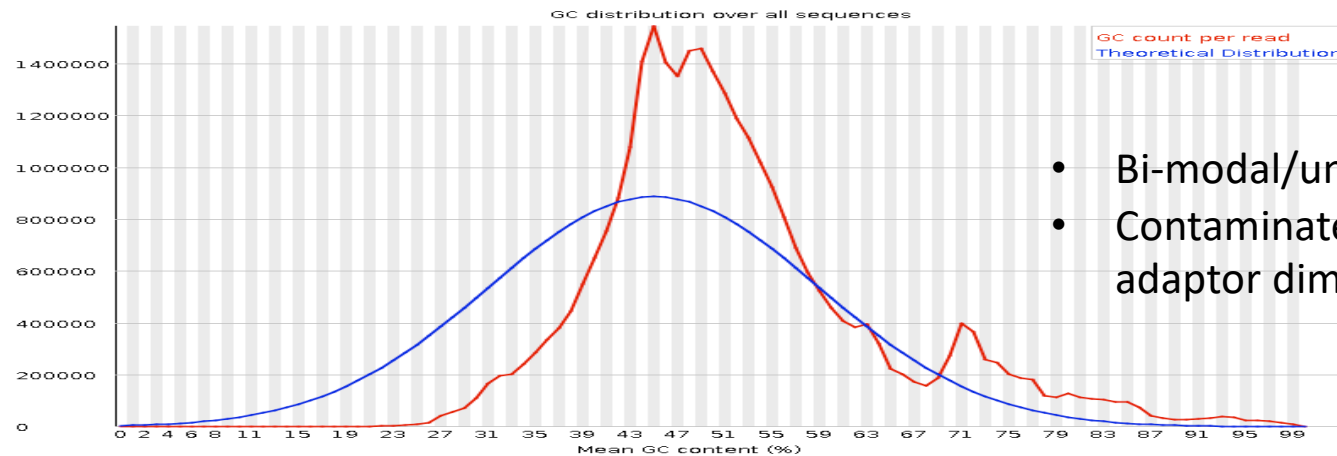


- Bi-modal/unusual distribution
- Contaminated/biased subset, i.e. adaptor dimmers, rRNA etc

Sequence duplication - FastQC

- Relative number of sequences with different degrees of duplication

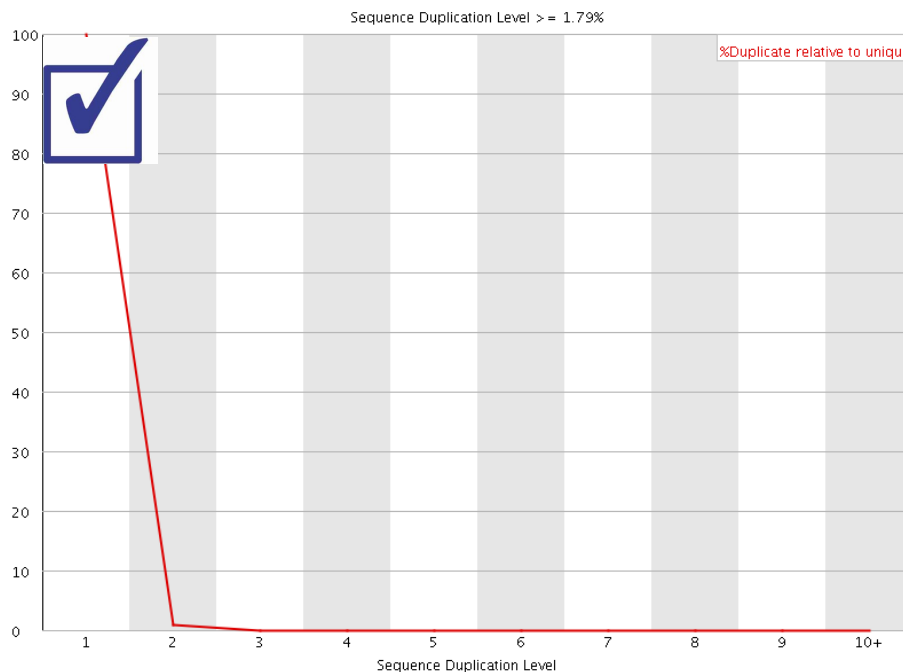

- Essentially no duplication

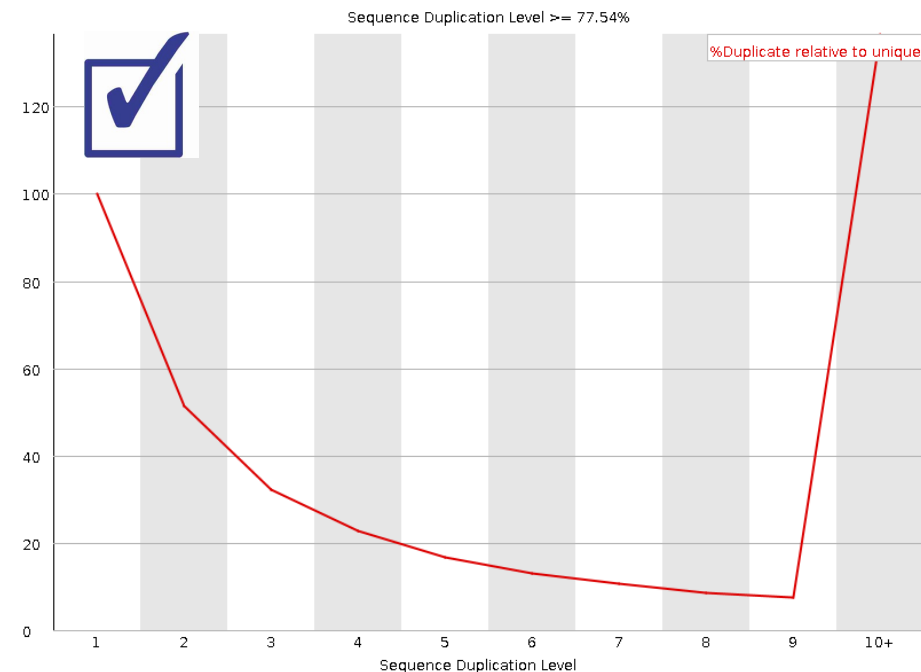

High duplication levels:

- DNA-seq: PCR over amplification, too little input material
- Normal in RNA-seq: high expression

Overrepresented sequences - FastQC



- Sequences make up >0.1 % of the total
- Compare those with a contamination database for finding contamination (i.e. adaptor dimmers)

Overrepresented sequences

Sequence	Count	Percentage	Possible Source
GGAAGAGCACACGTCTGAACTCCAGTCACCAGATCATCTCGTATGCCGTC	75874	1.5613887498682963	TruSeq Adapter, Index 7 (100% over 50bp)
GGAAGAGCACACGTCTGAACTCCAGTCACCGATGTATCTCGTATGCCGTC	7636	0.15713900010536297	TruSeq Adapter, Index 2 (100% over 50bp)
GGAAGAGCACACGTCTGAACTCCAGTCACACAGTGATCTCGTATGCCGTC	7539	0.1551428656095248	TruSeq Adapter, Index 5 (100% over 50bp)
GGAAGAGCACACGTCTGAACTCCAGTCACGCCAATATCTCGTATGCCGTC	5117	0.10530123933199874	TruSeq Adapter, Index 6 (100% over 50bp)


- Can be normal and biologically meaningful
 - highly expressed transcripts
 - high copy number repeats
 - Less diverse library (amplicons)

Adapter Content - FastQC

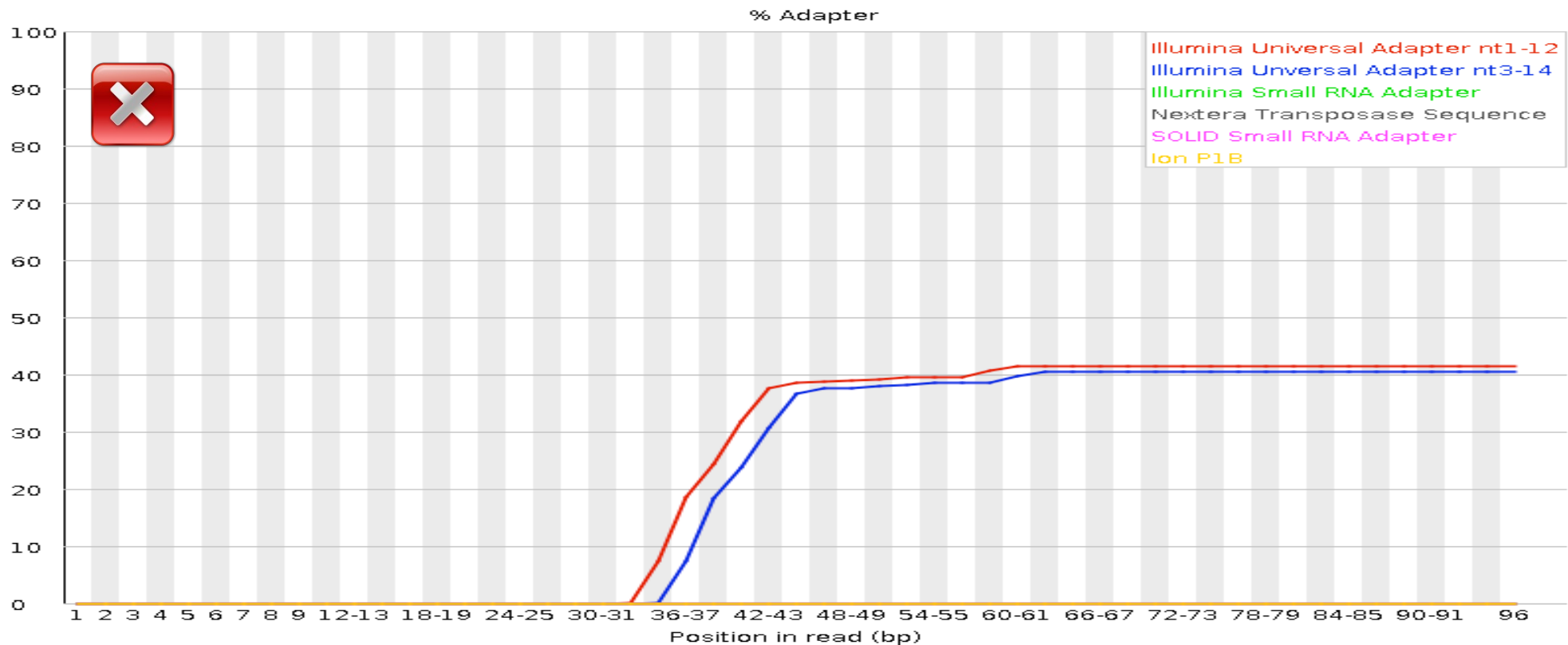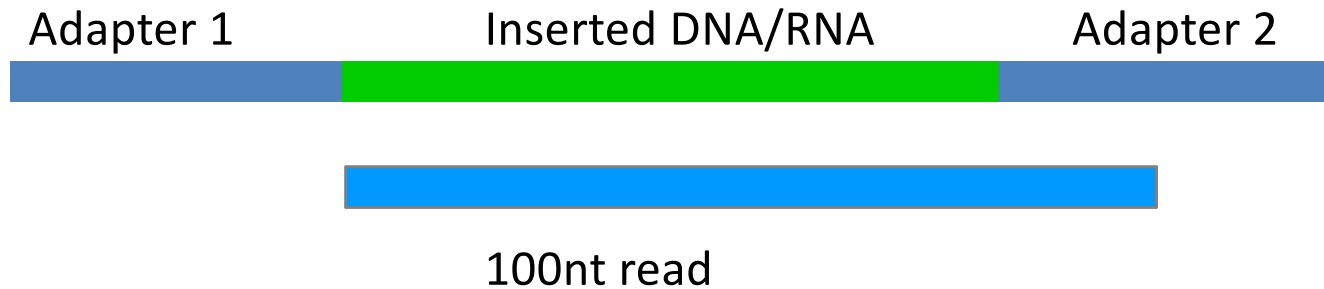


- How accurate was the sequencing → Fastqc
- Are these reads the intended ones → FastqScreen

Contamination Check - FastqScreen

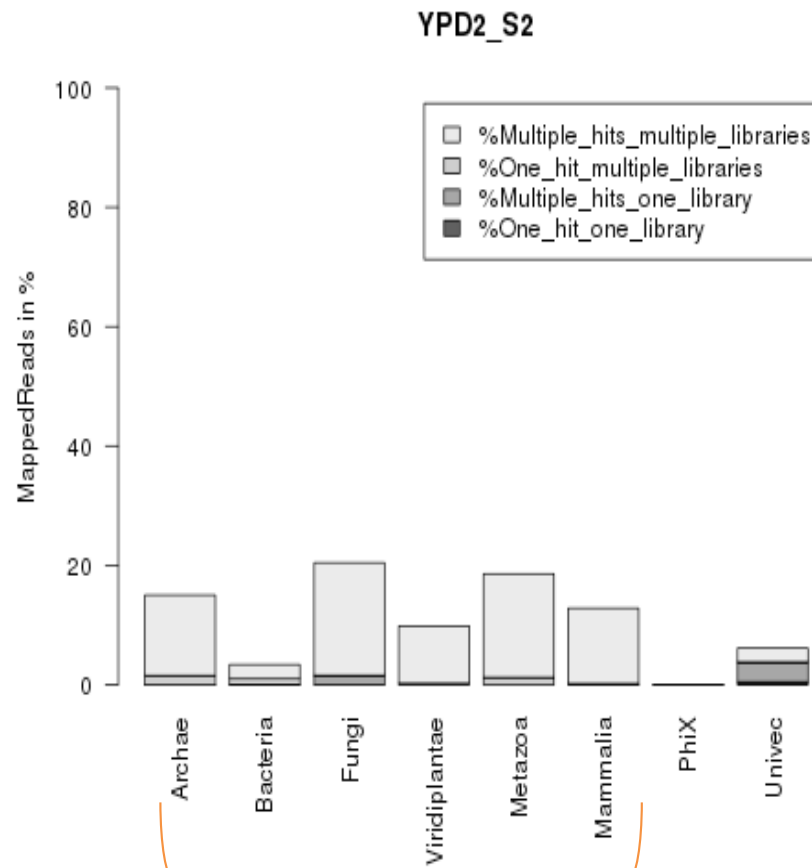

rRNA-Content (Silva)

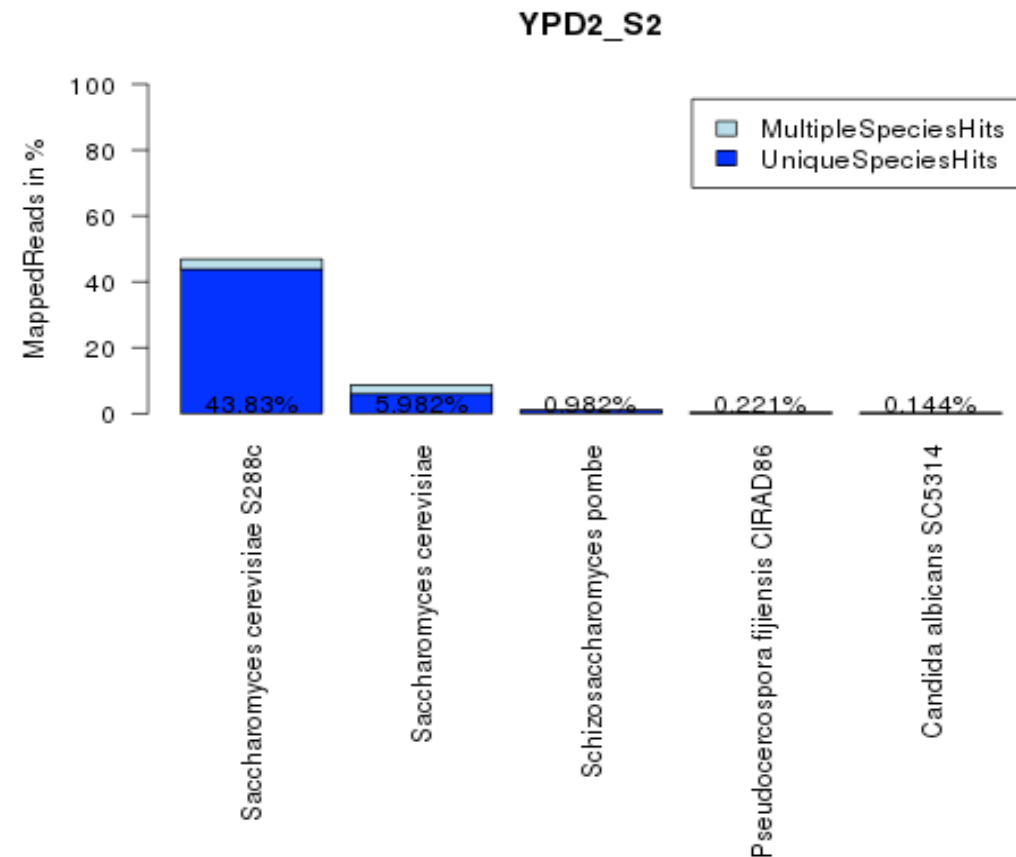

Mapping to RefSeq mRNAs
(all species)

Data preprocessing common tasks

1. Trimming: remove bad bases from (end(s) of) reads

- Adapter sequence
- Low quality bases

2. Filtering: remove bad reads

- Low quality reads
- Contaminating sequences
- Low complexity reads (repeats)
- Short (<20bp) reads – they slow down mapping software

Data preprocessing software

- fastp
 - <https://github.com/OpenGenome/fastp>
 - Adapter trimming, quality trimming & filtering, ...
- Trimmomatic
 - <https://github.com/usadellab/Trimmomatic>
 - Adapter trimming, quality trimming & filtering, ...
- FlexBar (FAR)
 - <https://github.com/seqan/flexbar>
 - Flexible barcode detection and adapter removal
- FASTX
 - http://hannonlab.cshl.edu/fastx_toolkit/
 - Reformat, stats, collapse duplicated reads, trim, filter, reverse compliment
- TagCleaner
 - <http://tagcleaner.sourceforge.net>
 - Trim MIDs or adaptors, demultiplexing
- DeconSeq
 - <http://deconseq.sourceforge.net>
 - Remove potential contaminants

Recommendations

- Always generate quality control plots visualizing key characteristics for all libraries
- Trim and/or filter data if needed
- Applications where erroneous reads are of concern:
 - de novo assembly
 - low coverage variant calling
- Applications that are more tolerant to low quality bases
 - RNA-seq